



InSite™ | *Listeria* | Rapid Environmental *Listeria* Species Test

Hygiene™ InSite *Listeria* is a quick and easy test for the detection of *Listeria* species in food processing environments. The test contains a liquid media formulated with growth enhancers, selective agents, and chromogenic compounds specific to the growth of *Listeria* species.


Benefits:

- Easy to use with no preparation steps required
- Large sponge tip maximizes sample collection for accurate results
- Eliminates costly lab fees for an equivalent test
- Color results are easy to interpret – brown/black indicates presence of *Listeria* species
- Presumptive positive results in as little as 24 hours
- Negative results in 48 hours
- Self-contained tube eliminates risk of contamination
- 12-month shelf life at refrigerated temperature (2 – 8 °C)

Sensitivity:

- <10 CFU *Listeria spp.*
- Detects heat injured *L. monocytogenes* at low levels (10-50 CFU of broth), even in the presence of high numbers of competing organisms

Interpreting Results

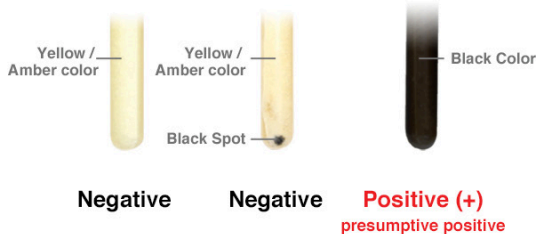

Catalog No.

Description

Quantity

IL050

InSite *Listeria* Environmental *Listeria spp.*

50

IL100


InSite *Listeria* Environmental *Listeria spp.*

100



Find support documents, instructional videos, and more at www.hygiene.com

InSite *Listeria* Procedure

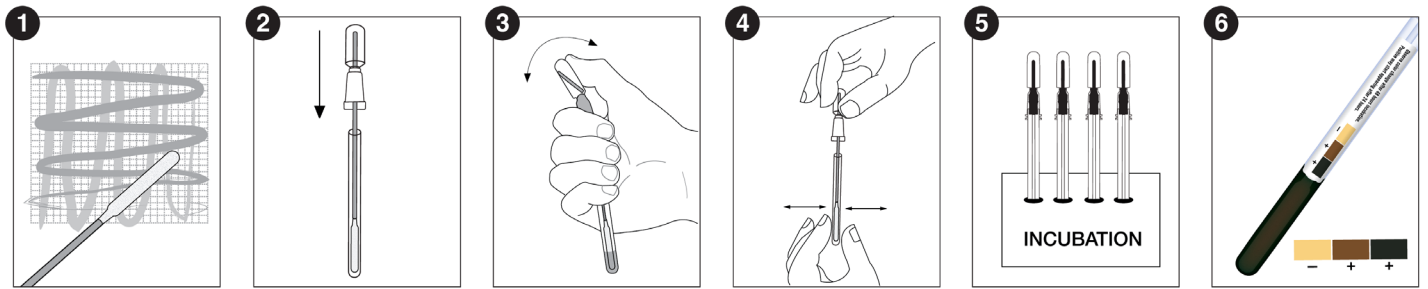

See an instructional demo at [YouTube HygieneTV](#) or [Twitter HygieneTV](#)

Related Products

Digital Dry Block Incubators

Digital incubators are capable of up to 105 °C and compatible with any Hygiene test device. Available in two sizes. Blocks are sold separately.




InSite *Salmonella*

InSite *Salmonella* is a quick and easy test for the detection of *Salmonella* species in 24 – 48 hours. The test contains a liquid media formulated with selective agents, growth enhancers, and color changing compounds specific to the growth of *Salmonella* species.



MicroSnap™ Coliform & *E. coli*

MicroSnap™ Coliform & *E. coli* are rapid tests for detection and enumeration of Coliform bacteria. Single figure organisms can be detected in the same day.

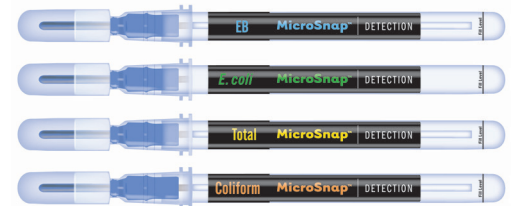

MicroSnap™ Total

MicroSnap Total delivers total viable count (TVC) results in the same day.

BAX® System Instruments

The BAX® System harnesses the power of DNA-based foodborne pathogen detection to help you make product decisions with confidence. An extensive assay portfolio, multiple platforms, and flexible protocols enable food labs of any size to take advantage of this benchmark technology.



Catalog No.	Description	Quantity
INCUBATOR/INCUBATOR2	Digital Dry Block Incubators	1
IS050	InSite <i>Salmonella</i> Environmental <i>Salmonella</i> Species Test	50
MS2-COLIFORM/MS2-ECOLI	MicroSnap™ Coliform & <i>E. coli</i>	100
MS2-TOTAL	MicroSnap™ Total	100
MCH2007	BAX® System X5	1
ASY2018	BAX® System Q7	1



Find support documents, instructional videos, and more at www.hygiene.com